



# Fear conditioning

How EthoVision XT supports your behavioral research



A white paper by Noldus Information Technology

# FEAR CONDITIONING: A DEEPLY-ROOTED BEHAVIOR

The most used readout for fear conditioning is freezing.


Fear conditioning is a fundamental behavioral process seen across a variety of species, from rodents to humans. It's a classic form of Pavlovian learning where subjects associate a neutral stimulus (like a light or tone) with an aversive event (mild foot shock).

There are two widely used variants of fear conditioning:

- Cued fear conditioning, which links a stimulus with a shock;
- Contextual fear conditioning associates the environment with an unpleasant stimulus.

The most used readout for fear conditioning is freezing, which can be easily measured with the Ugo Basile Fear Conditioning System in combination with EthoVision XT <sup>[1-5]</sup>.

*Creating automated tests and reliable detection of freezing behavior are essential for reproducibility of your results. Noldus EthoVision XT combined with the Ugo Basile Fear Conditioning System provides a solution that is centered on just those two aspects.*

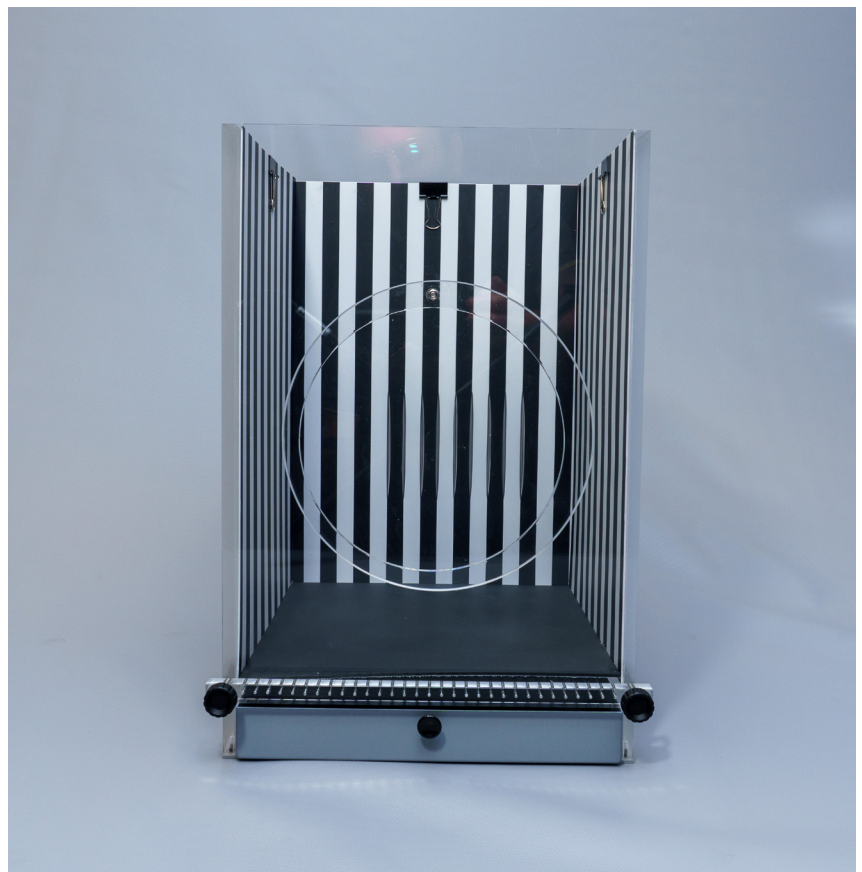

# HOW TO PERFORM A FEAR CONDITIONING TEST

To ensure precision and reproducibility, Noldus offers a fully integrated solution for measuring fear conditioning. EthoVision XT for automated behavioral tracking and the UGO Basile's Fear Conditioning Chamber together streamline every step of your experiment.


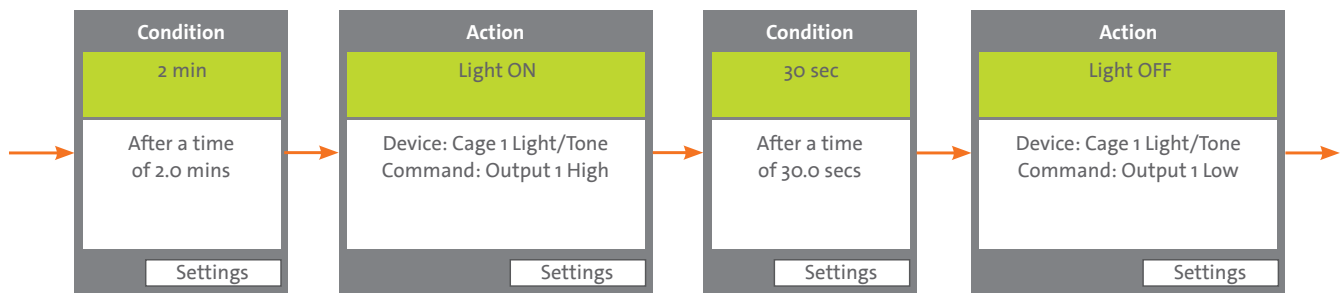


Figure 1. EthoVision XT offers an intuitive way to deliver stimuli for whichever protocol you decide to use.

Within the EthoVision XT software, you can start a protocol using a simple condition-action structure. The software automates your whole experiment with these simple steps. For example, after starting you can turn on the light for 30 seconds, after which a foot shock is administered. EthoVision XT correlates these triggers with the displayed behavior and gives you immediate insight into the behavior of your animals.

*EthoVision XT helps you creating routines to administer single or paired stimuli at specific times. You can also specify different protocols for different subjects tested simultaneously.*

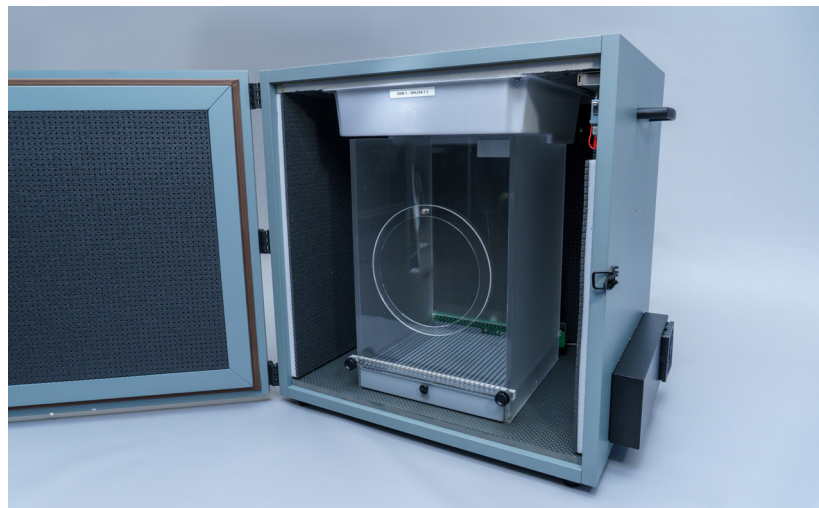


# DETECT FREEZING BEHAVIOR

EthoVision XT provides parameters that can accurately track mice for immobility that is nearly identical to scoring by human observation.

Freezing behavior is one of the most important readout parameters for fear conditioning research. EthoVision XT analyzes this behavior by tracking the amount of pixel changes and the change in shape of the animal. With frame-by-frame analysis the software can easily distinguish between an active state (while the animal is grooming, moving or exploring) and an inactive state. Pham *et al.*<sup>[6]</sup> used the parameter Mobility to score freezing automatically. According to the authors, EthoVision XT provides parameters that can accurately track mice and be used for automated scoring of immobility that is nearly identical to scoring by human observation.

However, with EthoVision XT you save a tremendous amount of time and therefore money conducting your experiments. Furthermore, with automated behavioral analysis your experiment is not biased towards the individual observer.

Figure 2. The Activity parameter in EthoVision XT: Purple pixels indicate the subject is active (a) or immobile (b).

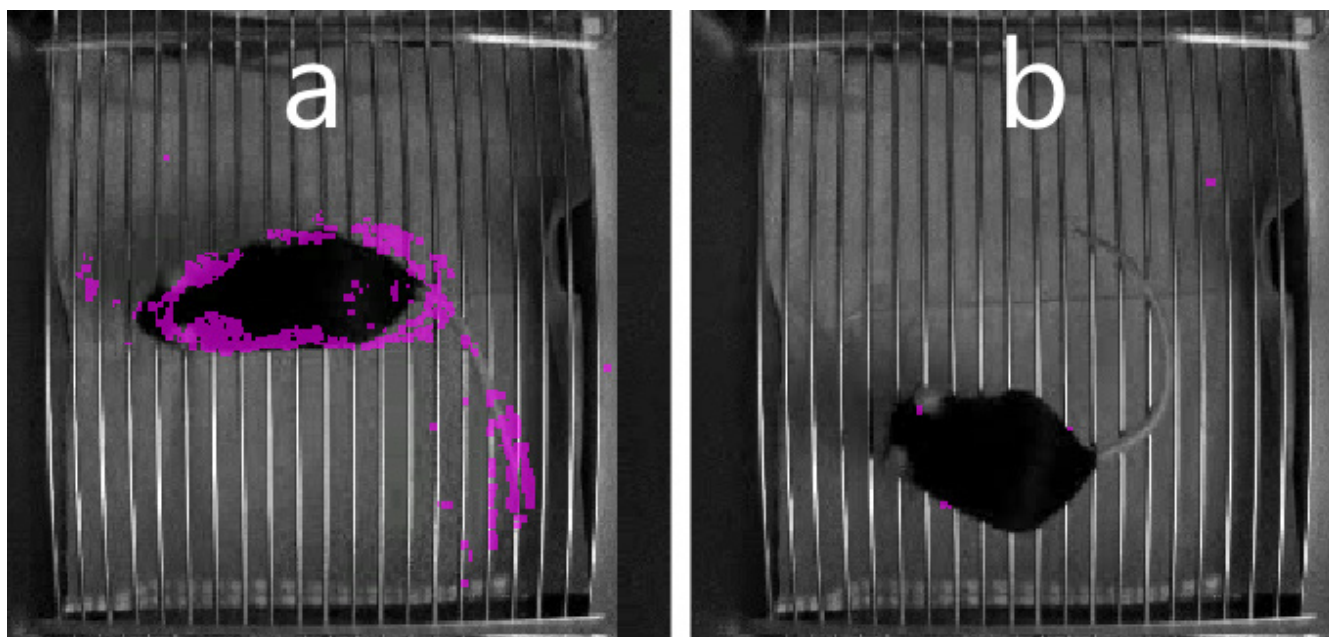

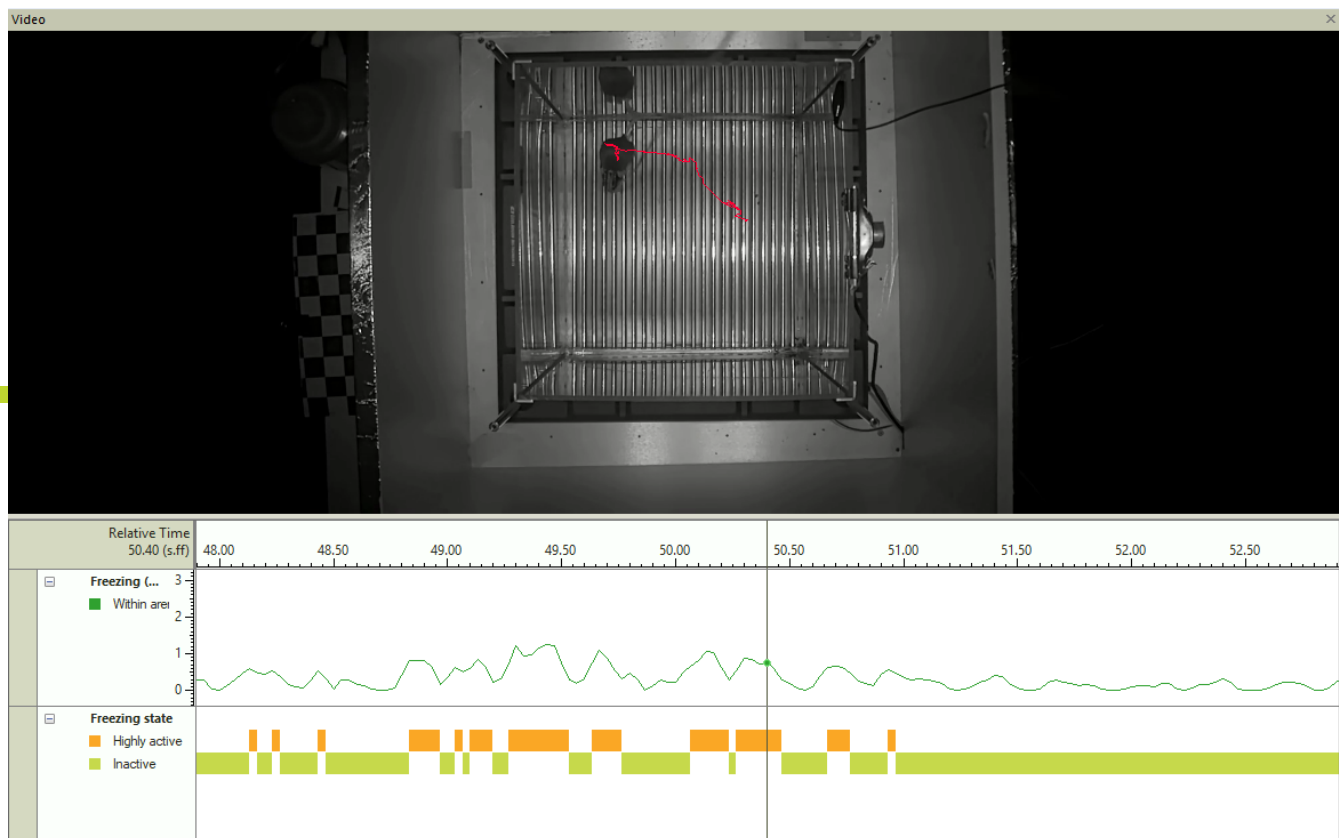


Figure 3. Activity state over time. Periods of inactive behavior could be classified as freezing.



# FEAR CONDITIONING IN PRACTICE

EthoVision XT is the best validated system for measuring fear conditioning.

## EFFECTS OF GENE OVEREXPRESSION ON FEAR LEARNING

Glutaminase C is an enzyme involved in the process that converts glutamine to glutamate, the most important neurotransmitter in excitatory synapses. A number of neurological diseases have been associated with changes in expression of the Glutaminase C gene. Wang *et al.*<sup>[8]</sup> generated transgenic mice that overexpressed Glutaminase C in the brain and tested their ability to learn fear. Compared with control mice, transgenic mice Nes-GAC mice had a significantly lower percentage of freezing compared on the last day of Cued fear conditioning, providing evidence that overexpression of Glutaminase C in the brain has deleterious effects on learning.

## BLOCKING INTERLEUKIN 1 SIGNALING PREVENTS STRESS-ENHANCED FEAR LEARNING.

There is a general consensus that post-traumatic stress disorder involves substantial immune system dysregulation. Exposure to severe stress is associated with an increase of the signaling molecule Interleukin- $1\beta$  in the hippocampus. Meghan E. Jones and collaborators (Jones *et al.*<sup>[9]</sup>) asked whether by simply blocking the interleukin signaling was sufficient to prevent fear learning. The authors used an antagonist to block Interleukin-  $1\beta$  in the hippocampus. Animals that received severe stress in a specific context followed by the interleukin antagonist infusion showed significantly less freezing than animals that received the same stress followed by vehicle. Freezing was measured using the *activity* parameter in EthoVision XT. The authors concluded that infusion of interleukin antagonist prevented the expression of stress-enhanced fear learning.

## FEAR CONDITIONING ELICITS RAPID MICROSTRUCTURAL CHANGES IN THE BRAIN

There is substantial knowledge about the anatomical, functional and molecular pathways involved in fear conditioning. Establishment of persistent long-term fear memory requires structural changes in the synapses in specific regions of the brain, particularly the amygdala and hippocampal

projections. Abby Ding and collaborators (Ding *et al.*<sup>[10]</sup>) used in vivo magnetic resonance diffusor tension imaging (DTI) before and after fear conditioning. To record freezing behavior in the subjects, the authors used the EthoVision XT parameter Mobility which analyzes subtle changes in the contour of the animal viewed from the camera in the fear conditioning chamber. The DTI revealed rapid microstructural changes in the amygdala and hippocampus in vivo, as early as one hour after conditioning. The study confirmed that the two areas of the brain are the most important in the regulation of fear.

#### DEPLETION OF RHEB GTPASE SIGNALING MOLECULE ELICITS MEMORY DEFECTS IN MICE

Ras genes encode for signaling molecules that usually work as two-state switches. If activated, a Ras molecule triggers other molecules, and ultimately determine signal pathways within the cell. Rheb (Ras homolog enriched in brain) GTPase is one of such molecules; it interacts with cell targets that are known to be implicated in life-span extension, as well as in neurodegenerative diseases such as Alzheimer's disease. However, it is not clear how Rheb GTPase participates in aging and neurodegenerative disorders. Shahani and collaborators (Shahani *et al.*<sup>[11]</sup>) crossed mice with Rheb with transgenic mice expressing an inducible promoter which depletes Rheb molecules. The mutant mice were then subject to contextual fear conditioning and the occurrence of freezing was recorded. Mice showed the same amount of freezing as controls. However, they performed less in spatial memory paradigms compared to controls, suggesting that Rheb depletion elicited spatial memory deficits.

## We are here to help you out.

Keep reading to find out more or contact us now to discuss your research goals.

**Continue reading to learn more about:**

- ✓ How EthoVision XT solves your problems
- ✓ Integration with 3<sup>rd</sup> party solutions
- ✓ Experiences from satisfied users

LEARN MORE

CONTACT US



# REFERENCES

1. Wehner, J.M., Radcliffe, R.A. (2004). Cued and Contextual Fear Conditioning in Mice. *Current Protocols in Neuroscience*, 8.5C.1-8.5C.14. John Wiley & Sons, Inc.
2. LeDoux, J.E. (2000). Emotion circuits in the brain. *Annu. Rev. Neurosci.*, **23**, 155-184.
3. Blanchard, R.J., Blanchard, D.C. (1969). Crouching as an index of fear. *Journal of Comparative and Physiological Psychology*, **67**, 370-375.
4. Marren, S. (2008). Pavlovian fear conditioning as a behavioral assay for hippocampus and amygdala function: cautions and caveats. *European Journal of Neuroscience*, **28**, 1661-1666.
5. Fanselow, M.S., Sterlace, S.R. (2014). Pavlovian fear conditioning: function, cause, treatment. In: McSweeney, F.K., Murphy, E.S. (Eds.), *The Wiley Blackwell Handbook of Operant and Classical Conditioning*. John Wiley & Sons, Ltd., Oxford, 117-141.
6. Pham, J., Cabrera, S.M., Sanchis-Segura, C., Wood, M.A. (2009). Automated scoring of fear related behavior using EthoVision software. *Journal of Neuroscience Methods*, **178**(2), 323-326. doi:10.1016/j.jneumeth.2008.12.021.
7. Anagnostaras, S.G., Wood, S.C., Shuman, T., Cai, D.J., LeDuc, A.D., Zurn, K.R., Zurn, J.B., Sage, J.R., Herrera, G.M. (2010). Automated assessment of Pavlovian conditioned freezing and shock reactivity in mice using the VideoFreeze system. *Frontiers in Behavioral Neuroscience*, **4**, 158. doi:10.3389/fnbeh.2010.00158
8. Wang, Y., Li, Y., Zhao, R., Wu, B., Lanoha, B., Tong, Z., Peer, J., Liu, J., Xiong, H., Huang, Y., Zheng, J. (2017). Glutaminase C overexpression in the brain induces learning deficits, synaptic dysfunctions, and neuroinflammation in mice. *Brain, Behavior, and Immunity*, **66**, 135-145.
9. Jones, M.E., Lebonville, C.L., Panizza, J.E., Balentine, M.E., Reissner, K.J., Lysle, D.T. (2018). Hippocampal interleukin-1 mediates stress-enhanced fear learning: A potential role for astrocyte-derived interleukin-1 $\alpha$ . *Brain, Behavior, and Immunity*, **67**, 355-363.
10. Ding, A.Y., Li, Q., Zhou, I.Y., Ma, S.J., Tong, G., McAlonan, G.M., Wu, E.X. (2013). MR diffusion tensor imaging detects rapid microstructural changes in amygdala and hippocampus following fear conditioning in mice. *PLoS ONE*, **8**(1), e51704.
11. Shahani, N., Huang, W.-C., Varum, M., Page, D.T., Subramaniam, S. (2017). Forebrain depletion of Rheb GTPase elicits spatial memory deficits in mice. *Neurobiology of Aging*, **50**, 134-143.





**INTERNATIONAL HEADQUARTERS**

Noldus Information Technology bv  
Wageningen, The Netherlands  
Phone: +31-317-473300  
E-mail: [contact@noldus.com](mailto:contact@noldus.com)

**NORTH AMERICAN HEADQUARTERS**

Noldus Information Technology Inc.  
Leesburg, VA, USA  
Phone: +1-703-771-0440  
Toll-free: 1-800-355-9541  
E-mail: [info@noldus.com](mailto:info@noldus.com)

**REPRESENTATION**

We are also represented by a  
worldwide network of distributors  
and regional offices. Visit our website  
for contact information.

**NOLDUS.COM**